

**Comfort-2**

Flex

# THE THERAPY CHAIR

HIGHEST COMFORT, ATTRACTIVE DESIGN IN ADDITION TO FUNCTIONALITY AND SAFETY

- The Comfort-2 Therapy Chairs are designed to assist and facilitate the work of physicians and medical staff
- They increase the patient's comfort considerably and therefore help to relieve tension during treatment
- Several optional accessories for optimal choice



PROFESSIONAL SEATING SOLUTIONS BY

**DIGITERM**

MEDICAL EQUIPMENT

# COMFORT-2 FLEX THERAPY CHAIR

## WITH EXCHANGEABLE UPHOLSTERY

Multiple independent adjustment possibilities by 3+1 motors offer maximum patient comfort, plus a reduction in work for staff.

### Electronic adjustments

- backrest
- seat tilt
- legrest
- footrest (optional)

### Tilt-up armrest

The solid upholstered armrest with height adjustment makes it possible to achieve the most suitable and restful arm position for any physique. The multi-layered soft armrest cushion is completely sealed against liquids and ensures improved patient comfort.

Optionally the distance between the armrests can be adjusted to 48 or 60 cm.

### Exchangeable upholstery

Can be removed easily with a few movements for cleaning or maintenance. It has been designed according to the latest developments in orthopaedics and ergonomics.

### Motorised footrest

The patient can adjust the position of the upholstered footrest to the most suitable and relaxing position at any time.

### Entry position

Adjustability of the Integrated table and Tilt-up armrest allow the patient to easily enter or leave the chair



### Trendelenburg position

The backrest adjusting range ensures that the patient's head is lower than the circulatory centre point of the body so as to meet concerning Medical Standards



\* Equipped with Integrated multifunctional table

### Easy-to-use 11 buttons handset

Automatic buttons


- One button operation for "Entry", perfectly horizontal "Bed" and Trendelenburg position

- The head-down position can be reached easily from any previous one by using the red button.
- Especially designed for medical applications

- Long-life micro-switches are located under an exchangeable keypad foil
- IP 67 Waterproof housing



# AVAILABLE OPTIONAL AND COMPLEMENTARY ACCESSORIES

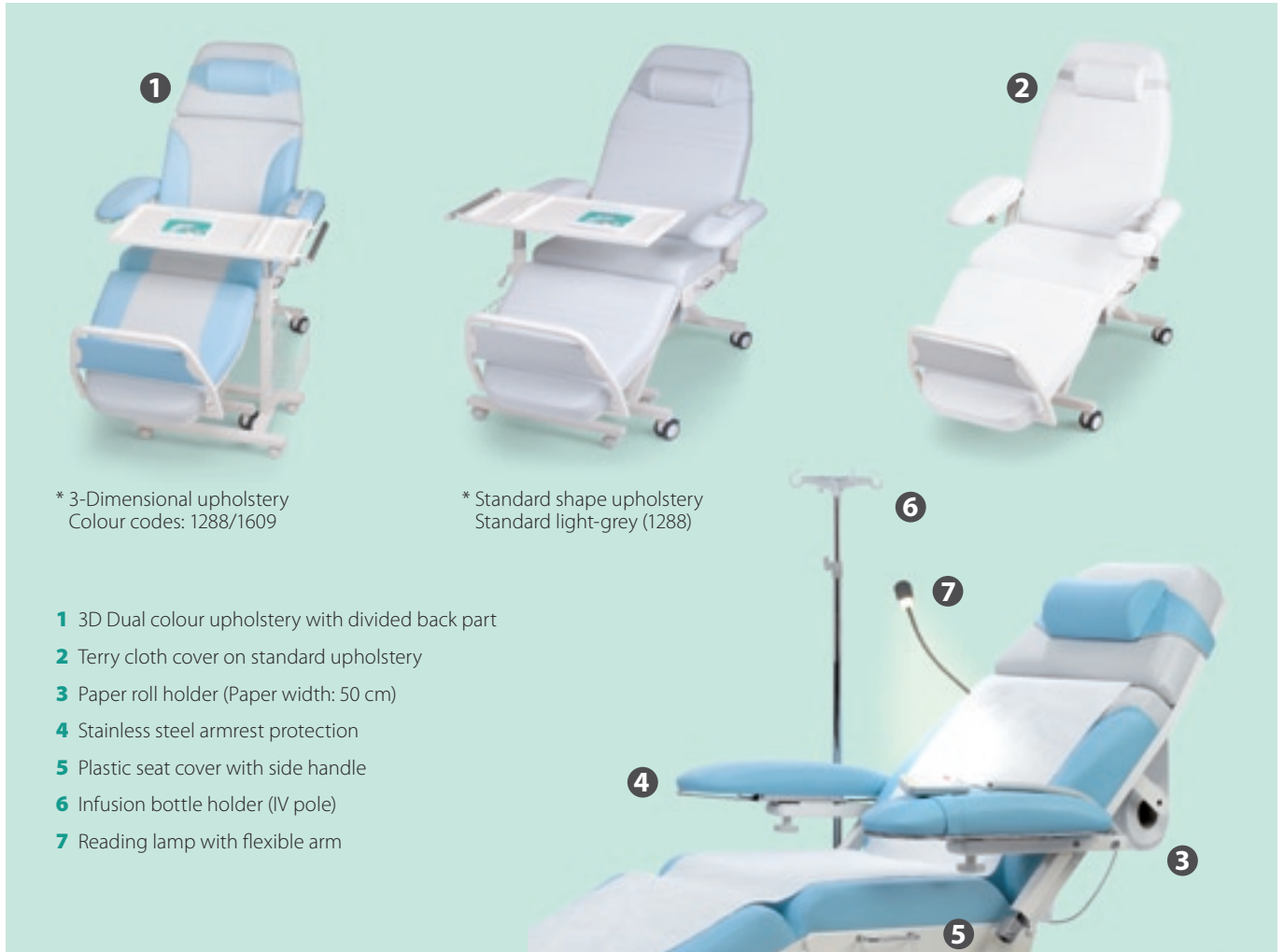

\* 3-Dimensional upholstery  
Colour codes: 1288/1609

\* Standard shape upholstery  
Standard light-grey (1288)

- 1 3D Dual colour upholstery with divided back part
- 2 Terry cloth cover on standard upholstery
- 3 Paper roll holder (Paper width: 50 cm)
- 4 Stainless steel armrest protection
- 5 Plastic seat cover with side handle
- 6 Infusion bottle holder (IV pole)
- 7 Reading lamp with flexible arm

## Comfort-T1 Bed-side table

Suitable for any type of COMFORT Therapy Chair and complement the latter during treatment.


The height of the table can be adjusted as required using the lever on its support stand. (75-98 cm) Size of trays: 40x56 cm and 40x17 cm.

## Integrated multifunctional table

provides space saving solution. It can be equipped on both sides of the chair.

The table can be positioned at an appropriate angle. The trays can be easily cleaned and lifted out.

## Sample Upholstery Colours

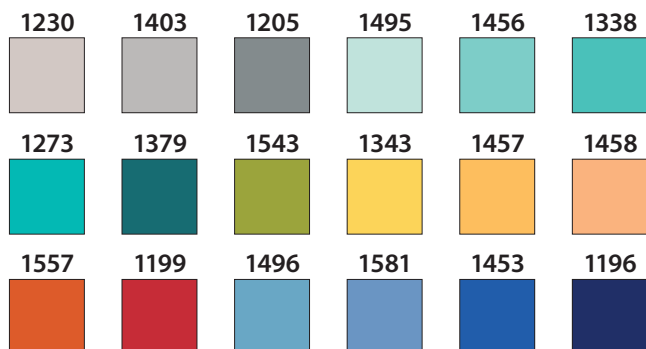

\* More than 30 colours are available upon request



Standard shape upholstery

# COMFORT-2 FLEX THERAPY CHAIR

ART.-NR. 7130382F

| Optional and complementary accessories                                                       | Article Nr. |
|----------------------------------------------------------------------------------------------|-------------|
| Comfort-T1 Bed-side table                                                                    | 6200501     |
| Integrated multifunctional table (built on the chair)                                        | 6200600     |
| Transparent plastic cover for the legrest                                                    | 7130240     |
| Transparent plastic cover for the footrest                                                   | 7130340     |
| Motorised footrest                                                                           | 7130404     |
| Reading lamp with flexible arm                                                               | 7130410     |
| Stainless steel armrest protection (It can be ordered only together with adjustable armrest) | 7130440     |
| Adjustable armrest (distance 48 or 60 cm)                                                    | 7130445     |
| Plastic seat cover with side handle                                                          | 7130460     |
| Terry cloth cover (standard white colour)                                                    | 7130501     |
| Infusion bottle holder (IV pole) with height adjustment                                      | 7130506     |
| Washable net bag for storing patient belongings                                              | 7130507     |
| Paper roll holder (It can be ordered only together with 3-Dimensional upholstery)            | 7130532     |
| Integrated table for Donor Scale                                                             | 7130580     |
| Standard shape upholstery in more than 30 colours (instead of standard light-grey)           | 7130690     |
| 3-Dimensional upholstery in single or dual-colour combination                                | 7130700     |

| General technical data          |                                                       |
|---------------------------------|-------------------------------------------------------|
| Backrest adjusting range        | 80° (-15° under horizontal in Trendelenburg position) |
| Seat tilt adjustment            | 15°                                                   |
| Legrest adjustment              | 35°                                                   |
| Footrest adjustment             | 29 cm                                                 |
| Degree of protection for motors | Protected against splashing water – IP66              |
| Armrest size (w x l)            | 17 x 57 cm                                            |
| Approved patient load           | max. 200 kg                                           |
| Size of castors                 | 4 x Ø7.5 cm with total brake                          |
| Size of upholstery              | 215 x 60 x 9,5 cm                                     |
| Weight of chair                 | 63 kg                                                 |
| Frame colour                    | RAL 9001                                              |
| Standard Power supply           | 230 V, 50/60 Hz (110 V, 127 V optional)               |
| Applied standards               | IEC 601; 93/42/EEC; CE                                |



**Digiterm Ltd.**  
Egri u. 30., H-3200 Gyongyos  
Tel.: +36 37 311 990  
Fax: +36 37 311 893  
E-mail: sales@digiterm.hu  
www.digiterm.hu

Approved by TÜV Product Service

