

**Comfort-4**

Eco

# THErapy CHAIR

HIGHEST COMFORT, ATTRACTIVE DESIGN IN ADDITION TO FUNCTIONALITY AND SAFETY

- The Comfort-4 Therapy Chairs are designed to assist and facilitate the work of physicians and medical staff
- They increase the patient's comfort considerably and therefore help to relieve tension during treatment
- Several optional accessories for optimal choice



PROFESSIONAL SEATING SOLUTIONS BY

**Digiterm**<sup>®</sup>  
medical equipment

\* Chair in the photo involves some optional elements

# COMFORT-4 ECO THERAPY CHAIR

## WITH VARIABLE HEIGHT

### Electronic adjustments

Multiple independent adjustment possibilities offer maximum patient comfort, plus a reduction in work for staff. (4 motors)

- backrest
- seat tilt
- legrest
- seat height

### Footrest adjustment

Stepless manually adjustable footrest provides the best comfort for any physique.

### Tilt-up armrest

The solid upholstered armrest makes it possible to achieve the most suitable and restful arm position for any physique.

It can be adjusted in height and entry is facilitated by simply tilting it up.

The distance between armrests can be adjusted also by swivelling completely. (48 or 60 cm)

The multilayered soft armrest cushion ensures improved patient comfort.

### Exchangeable upholstery

Can be removed easily with a few movements for cleaning or maintenance. It has been designed according to the latest developments in orthopaedics and ergonomics. Blood resistant artificial leather, which is resistant to disinfectants too.

### Entry position

Electronic height adjustment and Tilt-up armrest allowing the patient to easily enter or leave the chair



### Trendelenburg position

The backrest adjusting range ensures that the patient's head is lower than the circulatory centre point of the body in compliance with EN 60601 Medical Standard



\* Chairs in the photo involve some optional elements


### Easy-to-use 11 buttons handset

Automatic buttons

- One button operation for "Entry", perfectly horizontal "Bed" and "Trendelenburg" position

- The head-down position can be reached easily from any previous one by using the red button.
- Especially designed for medical applications

- Long-life micro-switches are located under an exchangeable keypad foil
- IP 67 Waterproof housing

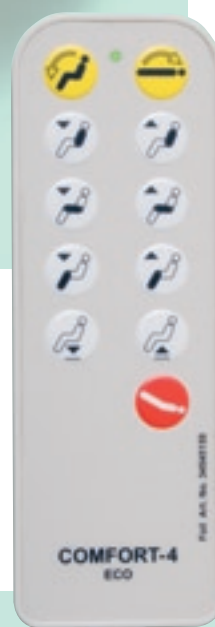

COMFORT-4  
ECO

REF. 845. 040. 30401110

# AVAILABLE OPTIONAL ACCESSORIES



3-Dimensional upholstery in dual-colour combination  
Colour codes: Perle 1288 / Dolphin 1609

Standard shape upholstery in standard light-grey colour  
Colour code: Perle 1288

- 1 3-Dimensional upholstery in dual-colour combination
- 2 Integrated multifunctional table
- 3 Central operating brake for quick fixing of 4 castors
- 4 Comfort-T1 Bed-side table
- 5 Terry cloth protecting cover on standard shape upholstery
- 6 Paper roll holder (Paper width: 50 cm)
- 7 Stainless steel armrest protection
- 8 Plastic side cover with handle on both sides of seat
- 9 Infusion bottle holder (IV pole) with height adjustment
- 10 Reading lamp with flexible arm

## Comfort-T1 Bed-side table

Suitable for any type of Comfort Therapy Chair and complement the latter during treatment.

The height of the table can be adjusted as required using the lever on its support stand. (75-98 cm) Size of trays: 40x56 cm and 40x17 cm.

## Integrated multifunctional table

provides space saving solution. It can be equipped on both sides of the chair.

The table can be positioned at an appropriate angle. The trays can be easily cleaned and lifted out.

## Sample Upholstery Colours

perle 1288	sand 1230	chrom 1205	mint 1495	aqua 1456	türkis 1338
agave 1273	limone 1543	limonade 1343	mais 1457	apricot 1458	karneol 1557
kirsche 1199	dolphin 1609	bleu 1496	sky 1581	atoll 1453	royal 1196

\* More than 50 colours are available upon request



# COMFORT-4 ECO THERAPY CHAIR

ART.-NR. 7130368E

Optional accessories	Art.-Nr.
Comfort-T1 Bed-side Table	6200501
Transparent protecting plastic cover for the legrest	7130240
Transparent protecting plastic cover for the footrest	7130340
Central operating brake with 4x diameter 10 cm twin wheel castors	7130402
Central operating total or directional lock with 4x diameter 12,5 cm twin wheel castors	7130413
Reading lamp with flexible arm	7130410
Stainless steel armrest protection (one pair)	7130440
Handle on both sides of seat (plastic side cover not included)	7130455
Plastic side cover with handle on both sides of seat (one pair)	7130460
Trendelenburg foot-switches on both sides of underframe	7130482
Terry cloth (Frottier) protecting cover for upholstered parts (standard white colour)	7130501
Infusion bottle holder (IV pole) with height adjustment	7130506
Washable net bag for storing Patient's belongings	7130507
Paper roll holder (recommended together with 3-Dimensional upholstery)	7130532
Paper roll 10 pcs (L 57m, W 50cm)	7130542
Integrated multifunctional table (built on the chair)	7130600
Individual colour for Standard shape upholstery, selectable from more than 30 different colours of TUNDRA	7130690
Special flame retardant high-tech artificial leather for Standard shape upholstery, selectable from more than 25 different colours of PALMA (meets M1 and some other Standards) no additional colour surcharge	7130691
Special flame retardant high-tech artificial leather for 3-Dimensional upholstery in single or dual-colour combination, selectable from more than 25 different colours of PALMA (meets M1 and some other Standards) no additional colour surcharge	7130693
Distributor's individual logo on top of backrest	7130695
3-Dimensional shape upholstery design that makes Standard shape upholstery even more comfortable and attractive. Available in single or dual-colour combination, selectable from more than 30 different colours of TUNDRA, no additional colour surcharge	7130700
Flexible height adjustable armrest for universal medical applications (one pair)	7130800
Integrated gas spring in the arm-support for smoother and safer operation (one pair)	7130810

## General technical data

Seat height adjustment	54 - 76 cm	Armrest size	L 57 x W 17 cm
Backrest adjusting range	80° (-13° under horizontal in Trendelenburg position)	Size of castors	4 x Ø 10 cm swivel castors with total lock
Seat tilt adjustment	20°	Weight of chair	78 kg
Legrest adjustment	30°	Maximal load	200 kg weight capacity
Footrest adjustment	29 cm (manually)	Standard Power Supply	230 V, 50/60 Hz (100 V, 120 V optional)
Dimensions in seat position	L 145 x W 92 or L 145 x W 80 cm*	Technical data of motors	24V DC, IP 66
Size of upholstery	L 215 x W 60 x H 9.5 cm	Applied directive and standards	93/42/EEC and harmonized standards (EN60601-1; EN60601-2-38; etc.)

\* Depend on armrest and footrest setting

**Digiterm**<sup>®</sup>  
medical equipment

**Digiterm Ltd.**

Egri u. 30., H-3200 Gyöngyös

Tel.: +36 37 311 990

Fax: +36 37 311 893

E-mail: sales@digiterm.hu

www.digiterm.hu

Approved by TÜV Product Safety

